

UCSF

UC San Francisco Previously Published Works

Title

Fluorescent amplification for next generation sequencing (FA-NGS) library preparation.

Permalink

<https://escholarship.org/uc/item/4x548860>

Journal

BMC genomics, 21(1)

ISSN

1471-2164

Authors

Chiniquy, Jennifer
Garber, Megan E
Mukhopadhyay, Aindrila
et al.

Publication Date

2020

DOI

10.1186/s12864-020-6481-8

Peer reviewed

METHODOLOGY ARTICLE

Open Access



Fluorescent amplification for next generation sequencing (FA-NGS) library preparation

Jennifer Chiniquy^{1,2†}, Megan E. Garber^{1,3†}, Aindrila Mukhopadhyay¹ and Nathan J. Hillson^{1,2*} 

Abstract

Background: Next generation sequencing (NGS) has become a universal practice in modern molecular biology. As the throughput of sequencing experiments increases, the preparation of conventional multiplexed libraries becomes more labor intensive. Conventional library preparation typically requires quality control (QC) testing for individual libraries such as amplification success evaluation and quantification, none of which occur until the end of the library preparation process.

Results: In this study, we address the need for a more streamlined high-throughput NGS workflow by tethering real-time quantitative PCR (qPCR) to conventional workflows to save time and implement single tube and single reagent QC. We modified two distinct library preparation workflows by replacing PCR and quantification with qPCR using SYBR Green I. qPCR enabled individual library quantification for pooling in a single tube without the need for additional reagents. Additionally, a melting curve analysis was implemented as an intermediate QC test to confirm successful amplification. Sequencing analysis showed comparable percent reads for each indexed library, demonstrating that pooling calculations based on qPCR allow for an even representation of sequencing reads. To aid the modified workflow, a software toolkit was developed and used to generate pooling instructions and analyze qPCR and melting curve data.

Conclusions: We successfully applied fluorescent amplification for next generation sequencing (FA-NGS) library preparation to both plasmids and bacterial genomes. As a result of using qPCR for quantification and proceeding directly to library pooling, the modified library preparation workflow has fewer overall steps. Therefore, we speculate that the FA-NGS workflow has less risk of user error. The melting curve analysis provides the necessary QC test to identify and troubleshoot library failures prior to sequencing. While this study demonstrates the value of FA-NGS for plasmid or gDNA libraries, we speculate that its versatility could lead to successful application across other library types.

Keywords: Library preparation, Next generation sequencing, NGS, SYBR green, Echo, High-throughput, qPCR

Background

Next generation sequencing (NGS) is becoming a pre-dominant tool in answering a broad range of biological questions. Its popularity can be attributed to its cost-effectiveness, its broad utility, and its multiplexing capabilities, which can be used to sequence hundreds, if not

thousands, of individual libraries simultaneously [1]. Because library preparation is step-intensive and cumbersome, multiplexing comes with a series of technical challenges [2]. Most notably, after DNA purification, each individual library requires individual quality controls (QC) to determine if the library amplified and the concentration after amplification. The concentration of each library is then used to determine the relative quantities of DNA so that they can be pooled in equal amounts. Accurate quantification of NGS libraries is important to ensure evenness of library pooling before sequencing. Current strategies available for quantifying

* Correspondence: njhillion@lbl.gov

[†]Jennifer Chiniquy and Megan E. Garber contributed equally to this work.

¹Biological Systems & Engineering Division, Lawrence Berkeley National Laboratory, Berkeley, CA 94720, USA

²U.S. Department of Energy Agile BioFoundry, Emeryville, CA, USA

Full list of author information is available at the end of the article



© The Author(s). 2020 **Open Access** This article is distributed under the terms of the Creative Commons Attribution 4.0 International License (<http://creativecommons.org/licenses/by/4.0/>), which permits unrestricted use, distribution, and reproduction in any medium, provided you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons license, and indicate if changes were made. The Creative Commons Public Domain Dedication waiver (<http://creativecommons.org/publicdomain/zero/1.0/>) applies to the data made available in this article, unless otherwise stated.

NGS libraries after amplification and purification include Quant-iT dsDNA Assay Kit (Thermo Fisher Scientific, Waltham, MA), real-time qPCR-based library quantification methods like the KAPA Library Quantification Kit (Roche, Pleasanton, CA) [3], and digital PCR [4–6]. Each of the aforementioned quantification methods have been validated for NGS library preparation, but are notably laborious when applying them to high-throughput workflows. Modifications to the standard library preparation workflow must be considered to make multiplexing a more viable option for researchers. Our goal was to establish a modified library preparation workflow that eases multiplexing by limiting the number of steps required per individual library. We accomplished this by bridging amplification and QC steps to bypass individual DNA purification using real-time quantitative PCR (qPCR). We call this library preparation workflow modification fluorescent amplification for NGS (FA-NGS).

A wide range of applications rely on qPCR, including genotyping analysis, medical diagnostics, gene expression profiling [7], phytopathogen identification [8], forensic studies [9], and validation of DNA microarrays [10]. In the context of NGS, qPCR has been used primarily for the quantification of NGS libraries, which is used to pool purified libraries in equal amounts, a crucial step in preventing poor quality sequencing data [11]. As opposed to employing quantification after DNA amplification and purification of individual libraries, we hypothesized that we could replace these two steps with a single qPCR step using SYBR Green I.

In addition to reducing the number of steps in the NGS library preparation workflow, qPCR with SYBR green I added supplementary benefits to the modified NGS workflow. qPCR is a widely used technique for nucleic acid detection and quantification that employs polymerases together with intercalating fluorescent dyes or optionally fluorescently labeled sequence-specific probes. The development of inhibitor-tolerant fusion polymerases such as the Pfu-Sso7d polymerase [12, 13], which is used with SYBR green I allows for flexibility in reaction conditions, including nucleic acid sample input types with NGS library preparation reagents. Other advantages include a wide dynamic range (up to 8 logs), and high sensitivity even with low volumes or low input quantities. Finally, because the data is collected in a closed-tube system, there is reduced risk for sample contamination [14].

A limitation of using PCR dyes such as SYBR is a lack of binding specificity. These nonspecific dyes, unlike probe-based assays, intercalate with any dsDNA including primer-dimer, which can lead to false positives. However, since the melting temperature of primer dimer is typically much lower than the amplicons of interest, melting curve analysis (MCA) at the end of qPCR enables easy detection of amplicon over primer dimers. We hypothesized that we could use MCA, enabled by


qPCR, to determine if individual libraries were properly amplified without any additional reagents or costs.

To evaluate whether a combination of qPCR and MCA could be applied to NGS library construction, we modified two distinct library workflows. We first tested FA-NGS with Illumina's Nextera XT (Illumina, San Diego, CA). The standard workflow uses enzymatic fragmentation with transposases followed by conventional PCR amplification with indexed primers, purification, and quantification. Our workflow modifications included replacing PCR and the final quantification with low volume qPCR and SYBR Green I, and adding MCA. We also adopted additional modifications for the Nextera workflow, which have been previously shown to have utility for multiplexing NGS libraries such as automation using Echo acoustic liquid handling and reduction in reaction volumes [15–17]. We next evaluated FA-NGS with adapter ligation-based library construction kit using the NEBNext Ultra II DNA library preparation kit. The standard workflow includes DNA shearing, a multi-step process for adapter ligation, intermediate purification and quantification, and PCR with indexed primers followed by a final purification and quantification. As with the first NGS library workflow, our modifications included replacing PCR and the final quantification with low volume qPCR and SYBR Green I, and adding MCA.


The libraries from both modified workflows utilizing FA-NGS (Fig. 1, Additional file 1: Fig. S1) were sequenced in this study. These modifications reduced the hands-on time to construct the library, the PCR reaction volume, and the overall risk of sample contamination [14]. We found that by replacing conventional PCR with qPCR in NGS library preparation, we were able to successfully use quantified measurements of fluorescence as a proxy for relative concentration to pool 96 individual libraries (Fig. 1, Additional file 1: Fig. S1). Because we replaced PCR with qPCR, we were also able to perform MCA as an intermediate QC to confirm library amplification. This was done with a simple procedure and did not incur additional reagent cost.

Results

We chose to apply FA-NGS to two commonly used methods for preparing libraries. We used Nextera-based library preparations of plasmid DNA and adapter ligation-based library preparations of purified genomic DNA (gDNA). Using an input titration assay, we determined that the concentration of input adapter-ligated (AL) gDNA or plasmid could be tracked by FA-NGS with NGS primers (Fig. 2, Additional file 2: Fig. S2 and Additional file 3: Fig. S3). We observed that we could apply Nextera tagmentation and FA-NGS to plasmids with starting amounts as low as 1.5 pg per 7.5 μ L reaction. We also found that the detection limit of FA-NGS above a 1000 end relative fluorescence unit (RFU)

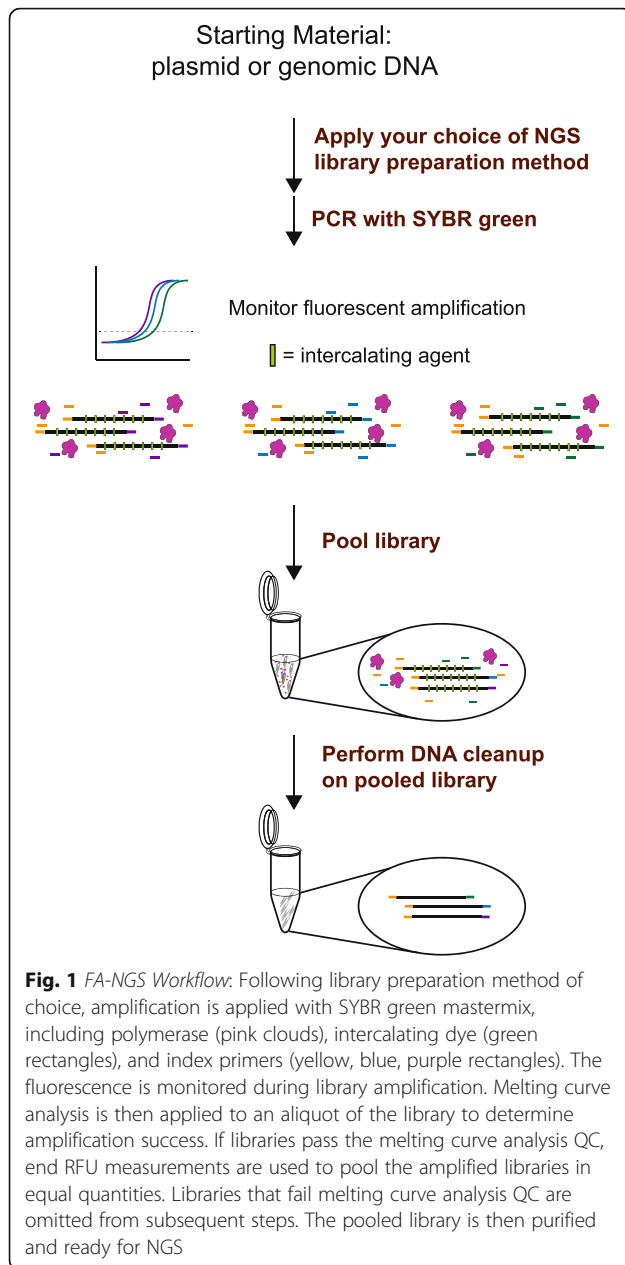

threshold for AL libraries was about 62.5 pg per 7.5 μ L reaction (Fig. 2a). For gDNA samples, the end RFU values scaled with the 2-fold dilution factor of the starting amount of the input (Fig. 2b). This indicated that the end RFU values could be used to estimate the transfer volumes required to generate a library of approximately equal concentrations of each indexed sample. MCA of the AL gDNA input showed that even samples with input below 62.5 pg still have amplified DNA, as observed by a high melting temperature peak (Fig. 2c, Additional file 4: Fig. S4B). For input values below 8 pg, no such high temperature peak was observed.

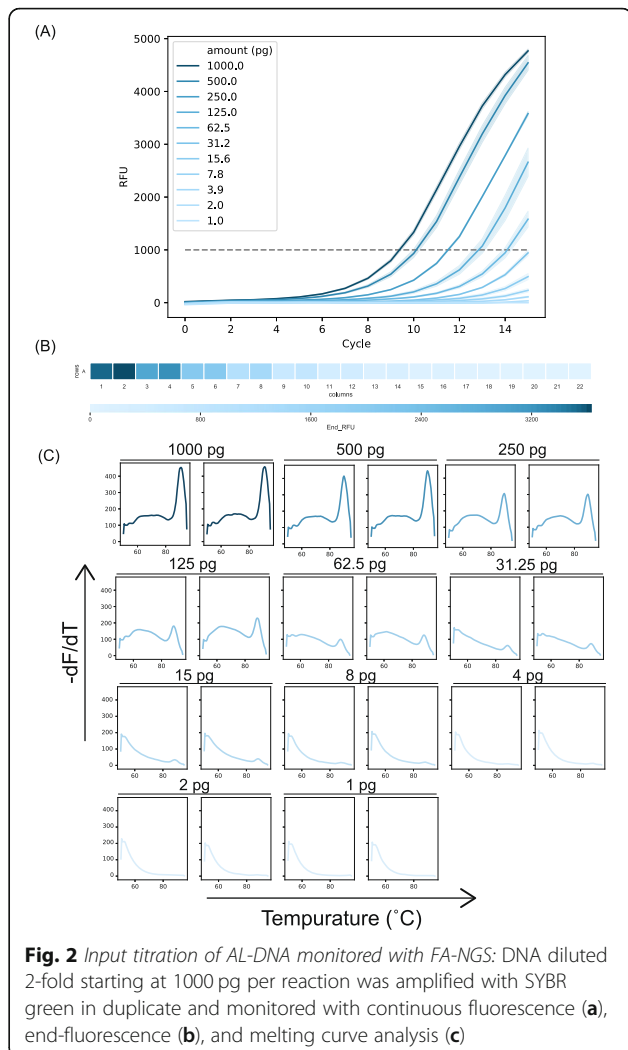



Combinations of dual index primers were used to uniquely index 96 samples. The plasmid library consisted of 4 plasmids (Table 1) with 24 replicates for each. The gDNA library was prepared with gDNA isolated from the same organism, *Pseudomonas stutzeri* RCH2, with 24 replicates of 4 starting amounts. The libraries were amplified and measured with continuous fluorescence. End RFU values were used as a proxy for relative concentration of each individual library. As expected, wells with similar starting amounts yielded similar end RFU values (Fig. 3b, Additional file 4: Fig. S4).

MCA (Additional file 5: Fig. S5 and Additional file 6: Fig. S6) was applied to assess the overall quality of each indexed sample. This enabled us to predict individual library failure prior to sequencing. As expected, a late melting peak was observed for wells with the appropriate template DNA. Samples without DNA template or with degraded DNA template were expected to have a low temperature peak, an indication of primer dimer and lack of library amplification (Fig. 2c, Additional file 3: Fig. S3).

Table 1 Plasmids used in Nextera library preparation. All plasmids used are available through the public instance of the ABF registry [18]. See the availability of data and materials section for additional information

Plasmid	Resistance	Size	Reference	Registry ID
pXMJ19	Chloramphenicol	6592	[19]	ABFPUB_000064
pGEN-292	Kanamycin	10,370	(Kirby et al. unpublished)	ABFPUB_000068
pskb3-CopR1598	Kanamycin	6029	[20]	ABFPUB_000072
pms6126	Carbenicillin	4038	[21]	ABFPUB_000070

To determine whether the end RFU values were a reasonable measurement of the final concentration of DNA after amplification, we sequenced the two libraries consisting of 96 combinations of dual indices with the Illumina MiSeq Reagent Nano kit. We observed that the percent of reads of each sample was comparable in both Nextera and AL libraries (Fig. 4). Both the Nextera and AL libraries yielded similar distribution of percent reads with a *p*-value of 1 (Additional file 7: Fig. S7). While each sample is not pooled at exactly 1.04% of the reads (expectation for optimal pooling from 96 samples), a majority of the samples from each library do fall below a 50% difference from optimal pooling range (Additional file 8: Fig. S8). Sequencing quality value scores for the PhiX Control Library and for the FA-NGS libraries were above the specification provided by Illumina of at least 80% Q30 [22] (Additional file 9: Fig. S9).

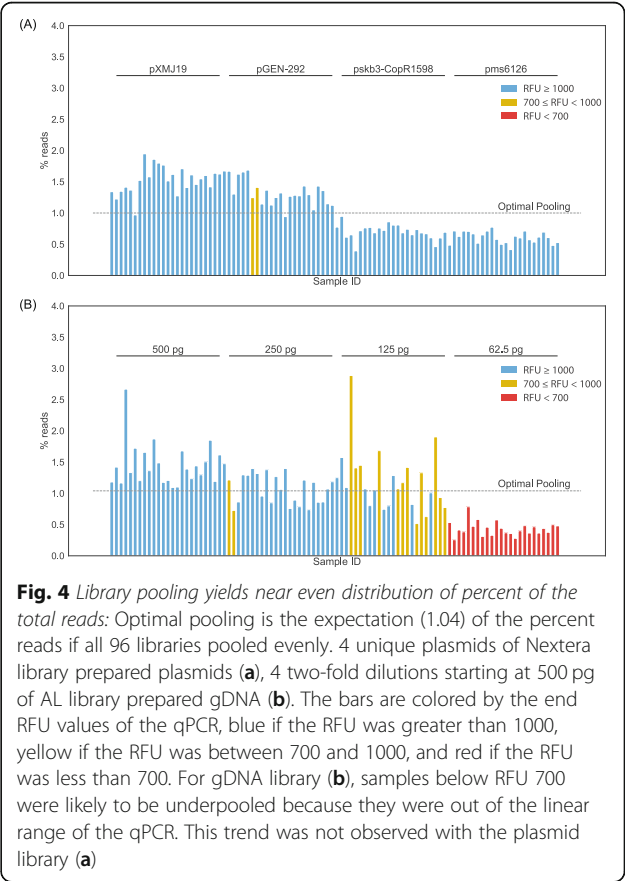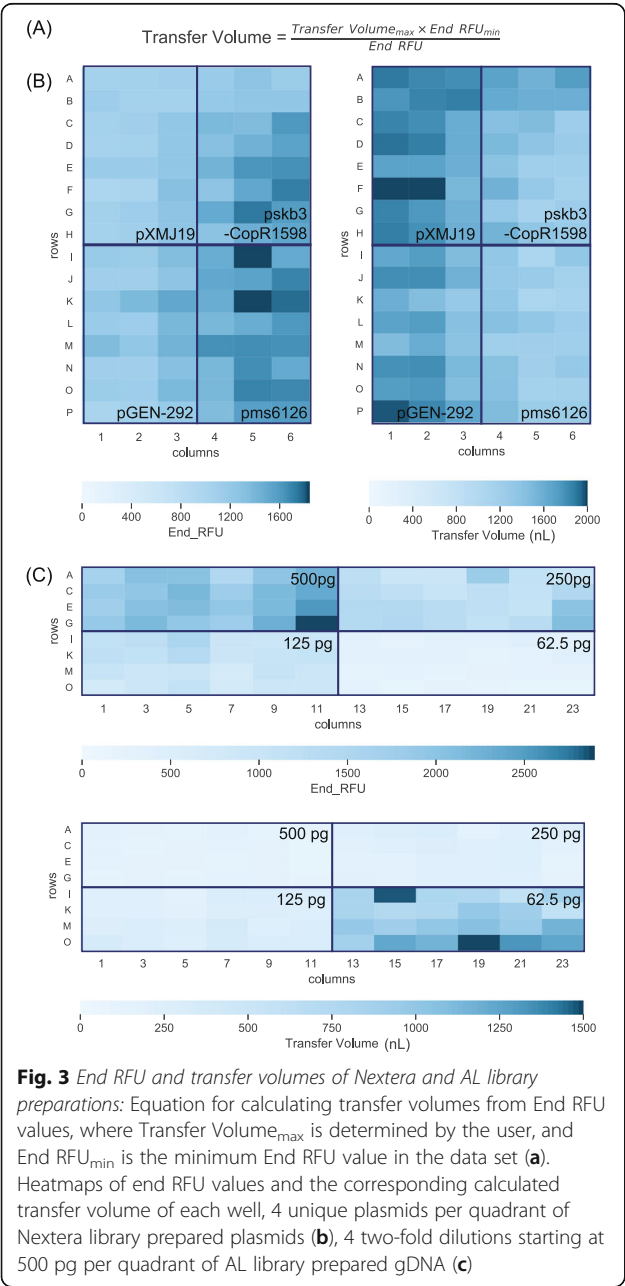

Discussion

In this report, we provide the validation of time-saving modifications to two common workflows for NGS library preparation. By implementing FA-NGS, we show that we can monitor amplification of each sample within an NGS library, perform quantification during qPCR to determine library normalization ratios for sequencing, and simultaneously perform a simple QC to confirm the amplified DNA template is present.

FA-NGS allows for accurate pooling directly after PCR, reducing the risk of user error. Typical NGS library preparation requires individual library purification prior to QC and pooling (Additional file 1: Fig. S1), which could lead to sample mix up or loss of library. Alternatively, unpurified and unquantified libraries could be pooled in equivolume amounts, however there would be no indication that the library amplified (increasing the potential for pooling primer dimer which would have deleterious effects during sequencing [23]) or that the samples varied in concentration. To avoid the latter issue, many PCR cycles (> 20 cycles) could be applied to fully saturate the reactions, however, this introduces the caveat of increased PCR bias, which can impact downstream interpretation of sequencing results. FA-NGS resolves these problems with a single step of amplification and quantification.

MCA enables users to perform individual library QC testing without DNA purification. Many workflows, especially Nextera-based library preparations lack QC tests until the final step of the protocol. When used on small aliquots of the amplified libraries, MCA confirms if an amplified library product is present. In applying FA-NGS, we suggest removing wells containing samples without amplified library, as determined by MCA, from the pooling calculations. Samples without amplified library may contain primer dimer, which if not removed in a final size selection can impact sequencing run quality and cause over clustering since smaller DNA fragments cluster more efficiently [24]. In accordance with our goal to make a user-friendly NGS library preparation workflow for multiplexed libraries, MCA adds an important QC step without additional costs or reagents.

The modified NGS workflow strategies described here are most beneficial for high-throughput NGS library applications, due to their compatibility with automation systems such as the Labcyte Echo. In order to make this modified workflow accessible to researchers with high-throughput NGS workflows, we made a FA-NGS software tool available on GitHub (see availability of data and materials) that can be used for set-up, analysis, and pooling of FA-NGS libraries. The “set-up tool” outputs a CSV file for primer transfer on the Echo acoustic liquid dispensing system. The “analysis tools” enable rapid determination of sample failure from end RFU and melting curve measurements. These tools output heatmaps

of the end-fluorescent values and small multiple MCA plots in the plate layout of the user’s choice. The “pooling calculator” outputs a CSV file, compatible with the Echo acoustic liquid dispensing system, for pooling libraries based on end RFU values.

A common goal in quantifying individual samples within a pooled library is to ensure that one given sample is not overrepresented as compared to another. Because of the potential variation in DNA shearing or fragmentation size distributions from NGS library preparation, it may be recommended to both measure the library concentrations and perform individual library size analysis to enable pooling in equimolar amounts. However, individual library size analysis for high-throughput sequencing applications can be laborious and costly. We demonstrate using relative fluorescence from qPCR to normalize library concentrations, without applying individual sizing analysis, that the distribution of sequencing reads (Fig. 4, Additional file 7: Fig. S7 and Additional file 8: Fig. S8) from using the FA-NGS library quantification strategy is sufficient to give comparable representation of 96 libraries.

To further validate the FA-NGS workflow, we performed sequencing quality assessment and included the Illumina PhiX Library Control for comparison. All libraries surpassed the recommended quality value score (Additional file 9: Fig. S9), per Illumina’s specification of the percentage of bases at or above Q30, indicating that the FA-NGS workflow successfully generated high quality sequencing reads.

While implementing FA-NGS has many benefits, the apparent pitfalls must be noted. We observed that samples in the gDNA library with end RFU values below 700 RFU were under-pooled (Figs. 3b and 4b). Because these samples were out of the fluorescence linear range of the qPCR, it is likely that the end-fluorescent values for these samples could not be used as a precise proxy for the relative DNA concentration and were therefore transferred at a volume that overestimated the relative concentration. We did not observe this trend with the plasmid library. We speculate that variation in percent reads for plasmid libraries could be an effect of amplifying beyond the exponential phase to the plateau phase of qPCR, such that the end-fluorescent values used for pooling were less precise (Additional file 10: Fig. S10). For standard library preparation workflows that use conventional PCR, there could be a comparable limitation during the final library quantification step, depending on the fluorescence linear range of the DNA quantification assay [25, 26]. For high-throughput applications employing FA-NGS, it may actually be preferable to allow amplification to the plateau phase, since initial sample concentrations or sample types may vary and the additional qPCR cycles help ensure all samples are sufficiently amplified.

The libraries in this report were sequence validated with Illumina Miseq Nano kits as proof of principle for our modified library workflow. Although the Nano kits

give relatively low read counts compared to other Illumina sequencing kits, they generated more than sufficient coverage for the 96 plasmids (Additional file 11: Table S1). However, for deep sequencing such as for bacterial genomes, larger MiSeq, HiSeq, or NovaSeq kits are generally recommended.

In this study, our modified workflow was successfully applied towards sequencing of full plasmids and bacterial genomes. These streamlined library preparations improve QC testing capability while reducing the total number of steps required to generate libraries.

Conclusions

Here we present a useful modification to conventional NGS library preparation workflows, FA-NGS, which was successfully incorporated into Illumina Nextera XT and NEBNext Ultra II DNA library preparation. We demonstrate the FA-NGS workflow ease of use with fewer overall steps than conventional library workflows, as well as an MCA QC test to confirm successful library construction before sequencing. An open source FA-NGS software tool is available to assist in implementing the workflow (see availability of data and materials). We expect that the application of FA-NGS will greatly benefit the production of any NGS library type which is amplified by PCR.

Methods

Plasmid DNA

Plasmid DNA was extracted from *E. coli* using the PureYield Plasmid Miniprep System (Promega, Madison, WI). *E. coli* cultures were grown overnight in LB medium supplemented with the appropriate resistance markers. DNA was quantified by the Qubit dsDNA HS Assay (Invitrogen, Carlsbad, CA). Each plasmid was diluted to 0.1 ng/μl for use in library construction. Plasmids are listed in Table 1. All plasmids used are available through the public instance of the ABF registry [18]. See the availability of data and materials section for additional information.

Enzymatic fragmentation and adapter ligation of plasmids

Plasmids were fragmented and ligated to partial Illumina adapters sequences in a single step using the Nextera XT DNA Library Preparation Kit (Illumina, San Diego, CA). Reactions were set up as described previously [27]. Using the Labcyte Echo 550 acoustic liquid dispensing system (Labcyte, Sunnyvale, CA), the Nextera fragmentation reaction volume was reduced to 1 μL total. Samples were incubated at 55 °C for 5 min for fragmentation and ligation, then at 70 °C for 15 min for heat inactivation of transposase enzymes. This reaction was immediately followed by qPCR.

Nextera library amplification of using real-time qPCR

Pairwise 8-nucleotide barcodes and the remainder of the Illumina adapter sequences were added using primers from the Nextera XT Index Kit v2. Primers were transferred using the Echo 550. Real-time qPCR and concurrent amplification and relative quantification of the library was facilitated on the CFX384 Touch Real-Time PCR Detection System using SsoAdvanced Universal SYBR Green Supermix (Bio-Rad, Hercules, CA). The Biomek Nx S8 was used to add SYBR Green and water. Reaction volumes were 7.5 μL total, with 0.25 μL each of the indexed Nextera primers, 3.75 μL of SYBR Green, 2.25 μL water, and 1 μL of the previous reaction. Cycling parameters were 72 °C for 3 min, 98 °C for 30 s, followed by 20 cycles of 10 s at 98 °C, 30 s at 63 °C and 3 min at 72 °C. The end RFU values were determined with CFX Manager software (Bio-Rad, Hercules, CA) for all libraries.

Shearing and adapter-ligation of genomic DNA

Pseudomonas stutzeri RCH2 was grown overnight in 5 mL UGA media [20]. UGA media contained 4.7 mM ammonium chloride, 1.3 mM potassium chloride, 2 mM magnesium sulfate, 0.1 mM calcium chloride, 0.3 mM sodium chloride, 5 mM sodium dihydrogen phosphate, 20 mM sodium lactate, and 25 mM MOPS. Vitamins and minerals were added as described by Widdel and Bak [28]. Genomic DNA was extracted using the Wizard genomic kit (Promega, Madison, WI) per the manufacturer's protocols. The resulting genomic DNA was diluted to 30 ng/μL and sheared with Covaris microTubes (Covaris, Woburn, MA) using the standard protocol for shearing DNA to 200 base pairs. Following shearing, 1 μg of sheared genomic DNA was ligated to adapters using the NEBNext Ultra II DNA library preparation kit (New England Biolabs, Ipswich, MA) according to the manufacturer's instructions without size selection and up until the final PCR. Size distribution of the adapter ligated DNA was verified on the Bioanalyzer (Agilent, Santa Clara, CA). It was then quantified with NanoDrop ND-1000 Spectrophotometer (NanoDrop Technologies, Inc., Wilmington, DE).

Adapter-ligation gDNA library amplification using real-time qPCR

7.5 μL PCR reactions were set up in a 384-well PCR plate (Bio-Rad, Hercules, CA) with 3.75 μL SYBR Green Supermix (Bio-Rad, Hercules, CA), 0.15 μL 50 μM i5 primer, 0.15 μL 50 μM i7 primer (Additional file 12: Table S2), and 3.45 μL (500, 250, 125 and 62.5 pg) adapter ligated DNA diluted in 1X Phosphate Buffer Saline adapter ligated DNA diluted in 1X Phosphate Buffer Saline. Cycling parameters were 72 °C for 3 min, 98 °C for 30 s, followed by 15 cycles of 10 s at 98 °C, 30 s at 65 °C and 1 min 45 s at 72 °C. The end RFU values were determined with CFX Manager software (Bio-Rad, Hercules, CA) for all libraries.

Set-up for melting curve and pooling for adapter-ligated libraries and Nextera libraries

The PCRs were diluted with 2.5 μ L water with the Bio-mek FX (Beckman Coulter, Indianapolis, IN). 7 μ L of the diluted PCRs were transferred to low dead volume (LDV) plates (Labcyte, Sunnyvale, CA) compatible with the Labcyte Echo 550 (Labcyte, Sunnyvale, CA) leaving 3 μ L of diluted PCR behind. The PCR plate was transferred to the CFX384 Touch Real-Time PCR Detection System (Bio-Rad, Hercules, CA) for MCA. The LDV plate was stored at -20°C until library pooling.

Melting curve analysis adapter-ligated DNA and Nextera libraries

The PCR plates were incubated in the CFX384 Touch Real-Time PCR Detection System (Bio-Rad, Hercules, CA) for MCA, gradually ramping from 50°C to 95°C using increments of 0.5°C , with 5 s at each temperature as fluorescence was monitored by CFX Manager software (Bio-Rad, Hercules, CA).

Library pooling calculation adapter-ligated DNA and Nextera libraries

The transfer volume of each PCR reaction was calculated with the end RFU values determined with CFX Manager Software (Bio-Rad, Hercules, CA). To calculate the transfer volume of a given well, the minimum fluorescence volume was multiplied by the maximum allowable transfer volume and was then divided by the actual fluorescence value of that well.

Library pooling

The LDV plate containing the diluted libraries was thawed to room temperature for library pooling. Each well was transferred to a 384-well PCR plate (Bio-Rad, Hercules, CA) with the Echo 550 (Labcyte, Sunnyvale, CA) using the transfer volumes calculated from the end RFU values. To avoid drip-back from the destination plate, the maximum volume transferred to each well in the 384-well PCR plate was 15 μ L. The contents of each transfer well in the destination plate were then pooled together in a 1.7 mL tube (Eppendorf, Hamburg, Germany). Following pooling, the library was cleaned up with AMPure beads (Beckman Coulter, Indianapolis, IN) according to the manufacturer's instructions. Library quality and size distribution was visualized on the Bioanalyzer (Agilent, Santa Clara, CA), and concentration was measured with Qubit dsDNA HS Assay Kit (Invitrogen, Carlsbad, CA).

Illumina sequencing

The Nextera and adapter-ligated libraries were sequenced with the MiSeq Reagent Nano kit v2 (Illumina, San Diego, CA), following Illumina's standard protocol. PhiX Control Library (v3) (Illumina, San Diego, CA) was

included with the libraries as an internal sequencing control. The Nextera library and adapter-ligated libraries ran for two rounds of 150 or 100 cycles, respectively.

Illumina sequencing analysis

Analysis of both sequencing runs was accomplished using embedded MiSeq Reporter (MSR) software (Illumina, San Diego, CA). The reads were aligned to the appropriate reference sequences with BWA-MEM [29]. Sequencing and alignment metrics were generated via MSR. For quality assessment of sequencing reads, BBTools Reformat [30] was used to generate quality value scores for individual libraries and the PhiX Control Library.

FA-NGS software tool

The FA-NGS software tool was written in python using pandas, numpy, seaborn, and matplotlib.pyplot libraries. The code is divided into four classes: set-up, plates, analysis, and pooling calculator. Set-up is used to output a CSV file with Labcyte echo instrument instructions for how to set up the multiplexed PCR plates with single or dual indexed primers. The input is an excel file with plate layouts for source plate, reverse primer destinations, and forward primer destinations. The plates class is used to customize plate setups. This software tool can use 96 well plates, 384 well plates, quadrants of a 384 well plate, and is suitable for some customizable setups. Analysis is used to read and visualize end RFU measurements and MCA directly from the Bio-Rad CFX instrument output files. The pooling calculator reads the end RFU measurements output file to determine pooling quantities. The calculator outputs instructions for pooling with Labcyte Echo instruments as well as a visualization of transfer volumes per well. The FA-NGS software tool was written in python and is available for download. See the availability of data and materials section for additional information.

Supplementary information

Supplementary information accompanies this paper at <https://doi.org/10.1186/s12864-020-6481-8>.

Additional file 1: Figure S1. Standard NGS and FA-NGS Workflow Comparison: Side-by-side comparison of the Standard NGS workflow (left) and modified FA-NGS workflow (right) highlights how FA-NGS can save time and hands on steps in preparing NGS libraries

Additional file 2: Figure S2. Input titration of Nextera library prepared plasmid monitored with FA: DNA diluted 4-fold starting at 100 pg per Nextera tagmentation reaction was amplified with SYBR green and monitored by continuous fluorescence (A), and melting curve analysis (B)

Additional file 3: Figure S3. MCA maxima correlates with input titration: Amount of input DNA is correlated with the local maxima of the MCA determined derivative RFU of the input titrations of both Nextera (A) and AL (B) library prepared samples, with R^2 equal to 0.807 and 0.842 respectively

Additional file 4: Figure S4. Distributions of end RFU and transfer volumes of AL libraries: End RFU values (A) and transfer volumes (B) of AL libraries generated with 500 (blue), 250 (red), 125 (yellow), and 62.5

(green) pg of input AL DNA are represented in a histogram. For adjacent histograms, * represents distributions with p -value < 0.001

Additional file 5: Figure S5. Melting curve analysis of Nextera prepared plasmids: The melting curve plot (temperature vs. negative derivative of fluorescence ($-dF/dT$)) of every well from the Nextera library is plotted. From left to right, the plasmids tested are pXMJ19, pskb3-CopR1598, pGEN-292, pms6126

Additional file 6: Figure S6. Melting curve analysis of AL prepared gDNA: The melting curve plot (temperature vs. negative derivative of fluorescence ($-dF/dT$)) of every well from the AL-gDNA library is plotted. From quadrant 1–4, the input concentrations are 500 pg, 250 pg, 125 pg, 62.5 pg

Additional file 7: Figure S7. Comparison of percent reads between Nextera and AL libraries shows similarities in output from two distinct NGS workflows: The distributions of Nextera (blue) and AL (red) libraries of percent reads are overlaid to highlight the similarities (p -value = 1) of sequencing output from these methods. The range of percent reads for the Nextera library (blue) was 0.39–1.95, with a mean of 1.04 and a standard deviation of 0.43. The range of percent reads for the AL library (red) was 0.25–2.89, with a mean of 1.04 and a standard deviation of 0.5

Additional file 8: Figure S8. Percent difference from sequence pooling of Nextera and AL libraries: The frequency of percent differences from the expected percent reads per sample (1.04) is represented as a histogram for the Nextera library (A), AL library (B)

Additional file 9: Figure S9. Sequencing quality scores of Nextera and AL libraries: The percentage of bases with \geq Q30 quality score for PhiX Control Library and for Nextera and AL libraries demonstrates sequencing quality for FA-NGS libraries

Additional file 10: Figure S10. Continuous fluorescence measurements of qPCR: RFU values per cycle number are plotted for 96 plasmid libraries (A) and 96 gDNA libraries (B)

Additional file 11: Table S1. Sequencing analysis for Nextera library: alignment analysis was performed using embedded MiSeq Reporter software

Additional file 12: Table S2. Primers for amplification of AL library

Abbreviations

FA-NGS: Fluorescent amplification for next generation sequencing; GDNA: Genomic DNA; MCA: Melting curve analysis; MSR: MiSeq Reporter; NGS: Next generation sequencing; QC: Quality control; QPCR: Quantitative PCR; RFU: Relative fluorescence unit

Acknowledgements

We would like to thank Cindi Hoover and Samuel Deutsch from the DOE Joint Genome Institute, who performed early development of the Nextera workflow. We would also like to thank Alexey Kazakov from Lawrence Berkeley National Laboratory, who assisted in assessment of sequencing quality for the adapter ligation workflow.

Authors' contributions

MG developed the software toolset for qPCR analysis and generation of robotic instructions. JC and MG both performed library preparation and sequencing. JC performed analysis of sequencing data from the MiSeq. AM and NJH supervised the work. All authors contributed to manuscript preparation and revision.

Funding

This work was part of the Agile BioFoundry (<http://agilebiofoundry.org>) supported by the U.S. Department of Energy, Energy Efficiency and Renewable Energy, Bioenergy Technologies Office, through contract DE-AC02-05CH11231 between Lawrence Berkeley National Laboratory and the U. S. Department of Energy. This work was also part of the ENIGMA, Ecosystems and Networks Integrated with Genes and Molecular Assemblies (<http://enigma.lbl.gov>), a Scientific Focus Area Program at Lawrence Berkeley National Laboratory and is supported by the U.S. Department of Energy, Office of Science, Office of Biological & Environmental Research under contract number DE-AC02-05CH11231 between Lawrence Berkeley National Laboratory and the U. S. Department of Energy. The funders had no role in study

design, data collection and interpretation, or the decision to submit the work for publication. The views and opinions of the authors expressed herein do not necessarily state or reflect those of the United States Government or any agency thereof. Neither the United States Government nor any agency thereof, nor any of their employees, makes any warranty, expressed or implied, or assumes any legal liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately owned rights. The United States Government retains and the publisher, by accepting the article for publication, acknowledges that the United States Government retains a nonexclusive, paid-up, irrevocable, worldwide license to publish or reproduce the published form of this manuscript, or allow others to do so, for United States Government purposes. The Department of Energy will provide public access to these results of federally sponsored research in accordance with the DOE Public Access Plan (<http://energy.gov/downloads/doe-public-access-plan>).

Availability of data and materials

The FA-NGS software tool was written in python and is available for download at: <https://github.com/AgileBioFoundry/FA-NGS>. All plasmids used are available through the public instance of the ABF registry: (<https://public-registry.agilebiofoundry.org/folders/2>) [18]. DNA sequencing was deposited in the Sequence Read Archive (SRA) database of the National Center for Biotechnology Information (NCBI) with Bioproject PRJNA599152.

Ethics approval and consent to participate

Not applicable.

Consent for publication

Not applicable.

Competing interests

NJH declares financial interests in TeselaGen Biotechnologies and Ansa Biotechnologies. All other authors declare that they have no competing interests.

Author details

¹Biological Systems & Engineering Division, Lawrence Berkeley National Laboratory, Berkeley, CA 94720, USA. ²U.S. Department of Energy Agile BioFoundry, Emeryville, CA, USA. ³Department of Comparative Biochemistry, University of California, Berkeley, CA 94720, USA.

Received: 22 October 2019 Accepted: 10 January 2020

Published online: 28 January 2020

References

- Park ST, Kim J. Trends in next-generation sequencing and a new era for whole genome sequencing. *Int Neurol J*. 2016;20(Suppl 2):S76–83.
- Neiman M, Sundling S, Grönberg H, Hall P, Czene K, Lindberg J, et al. Library preparation and multiplex capture for massive parallel sequencing applications made efficient and easy. *PLoS One*. 2012;7(11):e48616.
- Mardis E, McCombie WR. Library Quantification Using SYBR Green-Quantitative Polymerase Chain Reaction (qPCR). *Cold Spring Harb Protoc*. 2017(6):pdb.prot094714.
- Vogelstein B, Kinzler KW. Digital PCR. *Proc Natl Acad Sci U S A*. 1999;96(16):9236–41.
- Bhat S, Herrmann J, Armishaw P, Corbisier P, Emslie KR. Single molecule detection in nanofluidic digital array enables accurate measurement of DNA copy number. *Anal Bioanal Chem*. 2009;394(2):457–67.
- White RA, Blainey PC, Fan HC, Quake SR. Digital PCR provides sensitive and absolute calibration for high throughput sequencing. *BMC Genomics*. 2009; 10:116.
- Pastuszak-Lewandoska D, Bartosińska-Dyc A, Migdalska-Sek M, Czarnecka KH, Nawrot E, Domańska D, et al. HPV16 E6*II gene expression in intraepithelial cervical lesions as an indicator of neoplastic grade: a pilot study. *Med Oncol*. 2014;31(3):842.
- Schena L, Abdelfattah A, Mosca S, Li Destri Nicosia MG, Agosteo GE, Cacciola SO. Quantitative detection of *Colletotrichum godetiae* and *C. acutatum sensu stricto* in the phyllosphere and carposphere of olive during four phenological phases. *Eur J Plant Pathol*. 2017;149(2):337–47.

9. Goodwin C, Higgins D, Tobe SS, Austin J, Wotherspoon A, Gahan ME, et al. Singleplex quantitative real-time PCR for the assessment of human mitochondrial DNA quantity and quality. *Forensic Sci Med Pathol*. 2018; 14(1):70–5.
10. Ferrareso S, Vitulo N, Mininni AN, Romualdi C, Cardazzo B, Negrisolo E, et al. Development and validation of a gene expression oligo microarray for the gilthead sea bream (*Sparus aurata*). *BMC Genomics*. 2008;9:580.
11. Quail MA, Smith M, Coupland P, Otto TD, Harris SR, Connor TR, et al. A tale of three next generation sequencing platforms: comparison of ion torrent, Pacific Biosciences and Illumina MiSeq sequencers. *BMC Genomics*. 2012;13:341.
12. Wang Y, Prosen DE, Mei L, Sullivan JC, Finney M, Vander Horn PB. A novel strategy to engineer DNA polymerases for enhanced processivity and improved performance in vitro. *Nucleic Acids Res*. 2004;32(3):1197–207.
13. McInerney P, Adams P, Hadi MZ. Error rate comparison during polymerase chain reaction by DNA polymerase. *Mol Biol Int*. 2014;287430.
14. Kaltenboeck B, Wang C. Advances in real-time PCR: application to clinical laboratory diagnostics. *Adv Clin Chem*. 2005;40:219–59.
15. Mayday MY, Khan LM, Chow ED, Zinter MS, DeRisi JL. Miniaturization and optimization of 384-well compatible RNA sequencing library preparation. *PLoS One*. 2019;14(1):e0206194.
16. Minich JJ, Humphrey G, Benitez RAS, Sanders J, Swafford A, Allen EE, et al. High-Throughput Miniaturized 16S rRNA Amplicon Library Preparation Reduces Costs while Preserving Microbiome Integrity. *mSystems*. 2018; 3(6):00166–18.
17. Shapland EB, Holmes V, Reeves CD, Sorokin E, Durot M, Platt D, et al. Low-cost, high-throughput sequencing of DNA assemblies using a highly multiplexed Nextera process. *ACS Synth Biol*. 2015;4(7):860–6.
18. Ham TS, Dmytriv Z, Plahar H, Chen J, Hillson NJ, Keasling JD. Design, implementation and practice of JBEI-ICE: an open source biological part registry platform and tools. *Nucleic Acids Res*. 2012;40(18):e141.
19. Jakoby M, Ngouoto-Nkili C-E, Burkovski A. Construction and application of new *Corynebacterium glutamicum* vectors. *Biotechnol Tech*. 1999;13(6):437–41.
20. Garber ME, Rajeev L, Kazakov AE, Trinh J, Masuno D, Thompson MG, et al. Multiple signaling systems target a core set of transition metal homeostasis genes using similar binding motifs. *Mol Microbiol*. 2018;107(6):704–17.
21. Shih PM, Vuu K, Mansoori N, Ayad L, Louie KB, Bowen BP, et al. A robust gene-stacking method utilizing yeast assembly for plant synthetic biology. *Nat Commun*. 2016;7:13215.
22. Illumina, Inc. MiSeq system data sheet. 2018 [cited 2019 4];770–2011-001-B: 1–3. Available from: https://www.illumina.com/documents/products/datasheets/datasheet_miseq.pdf
23. Ravi RK, Walton K, Khosroheidari M. Miseq: a next generation sequencing platform for genomic analysis. *Methods Mol Biol*. 2018;1706:223–32.
24. Bourlat SJ, Hanel Q, Finnman J, Leray M. Preparation of amplicon libraries for Metabarcoding of marine eukaryotes using Illumina MiSeq: the dual-PCR method. *Methods Mol Biol*. 2016;1452:197–207.
25. Holden MJ, Haynes RJ, Rabb SA, Satija N, Yang K, Blasic JR. Factors affecting quantification of total DNA by UV spectroscopy and PicoGreen fluorescence. *J Agric Food Chem*. 2009;57(16):7221–6.
26. Haque KA, Pfeiffer RM, Beerman MB, Struewing JP, Chanock SJ, Bergen AW. Performance of high-throughput DNA quantification methods. *BMC Biotechnol*. 2003;3:20.
27. Thompson MG, Sedaghatian N, Barajas JF, Wehrs M, Bailey CB, Kaplan N, et al. Isolation and characterization of novel mutations in the pSC101 origin that increase copy number. *Sci Rep*. 2018;8(1):1590.
28. Widdel F, Bak F. Gram-Negative Mesophilic Sulfate-Reducing Bacteria. *The Prokaryotes*. Springer. 1992;3352–78.
29. Li H. Aligning sequence reads, clone sequences and assembly contigs with BWA-MEM. Oxford University Press. 2013;0(0):1–3.
30. Bushnell B. 2019. BBMap. Available from: <https://sourceforge.net/projects/bbmap/>

Publisher's Note

Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Ready to submit your research? Choose BMC and benefit from:

- fast, convenient online submission
- thorough peer review by experienced researchers in your field
- rapid publication on acceptance
- support for research data, including large and complex data types
- gold Open Access which fosters wider collaboration and increased citations
- maximum visibility for your research: over 100M website views per year

At BMC, research is always in progress.

Learn more biomedcentral.com/submissions

